

## Phoenix Image-Guided Laser System

### Compact, precise laser delivery

The Phoenix Research Labs' Image-Guided Laser System is user-friendly, producing precise, easy-to-deliver laser energy to facilitate research into the mechanisms of retinal damage and regeneration.

### Designed for small animals

The Image-Guided Laser System is coupled with the Micron IV Retinal Imaging Microscope, both of which are designed specifically for small animal research. No additional lab, bench space or clinical slit lamp delivery system is needed.

### Easy Targeting

Laser targeting using image-guidance is much easier and more precise due to the magnification of the bright field image. Simply align the animal to the Micron, position and focus the aiming beam, adjust the laser settings and deliver the energy.

### Clear advantages

With the cornea coupled to the objective lens using a gel interface, the eye is stabilized against movements from respirations and delivered power is more consistent than through open air. The coupling gel maintains hydration of the cornea, reducing the chances of media opacification.

### Publication quality

Use the bright field image as a guide to place the aiming beam on target areas. Easily capture photographs or videos for publications or presentations to strengthen research documentation. Integration with the Phoenix Micron IV Retinal Imaging Microscope enables the recording of pre/post laser treatment spots as well as documentation of longitudinal studies.

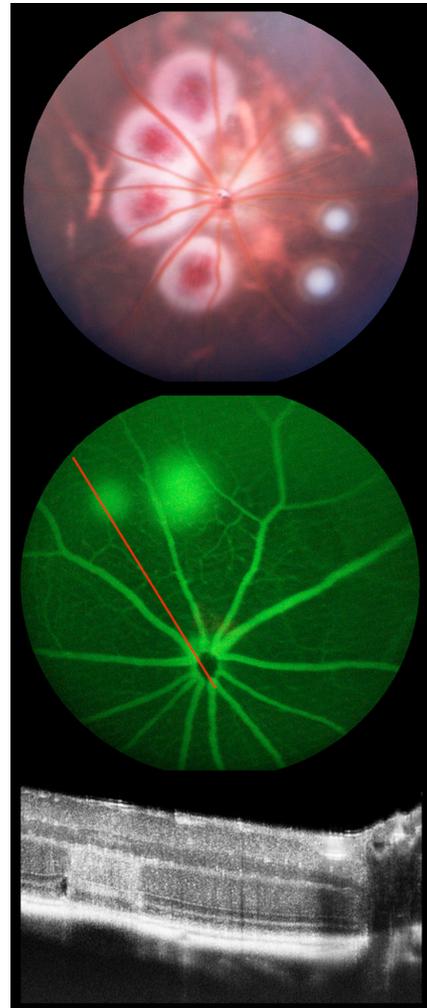


Image-Guided Laser System (top to bottom): Laser burns of varying intensities on a rat, fluorescein angiography on mouse with OCT beam overlay, corresponding OCT of penetration of laser burn on mouse.

Complete a comprehensive in-vivo retinal study at the spin of the filter wheel by performing fluorescein angiography to observe and document leakage patterns.

### Phoenix Image-Guided Laser fuels your research

The Micron IV design delivers a unique proposition for image-guided eye research tailored specifically for laboratory animals. The Phoenix Micron IV platform supports a family of additional turnkey research instruments including Slit Lamp Imaging, Focal ERG, and Image-Guided OCT. Unlike stand-alone adaptations of human instruments, the Phoenix range of products interface directly with the Micron IV to support your image-guided comprehensive studies. With a broad range of applications including basic research, toxicology, pharmaceutical efficacy testing and neurological research, the Micron IV is sure to fuel scientific discoveries, which is at the heart of the Phoenix mission.



The small, compact, easy to use Phoenix Image-Guided Laser makes research into retinal ablation simple.

#### SPECIFICATIONS

Laser wavelength	532 nm
Spot size	50 $\mu\text{m}$
Aiming beam laser wavelength	635 nm
Pulse duration	0-5 seconds
Maximum delivered laser power	0-600 mW
Laser Source	Meridian Merilas 532 green laser photocoagulator p/n 100470
Lenses	Separate standard Micron IV objective lenses for mouse and rat
	Separate laser objective lenses for mouse and rat
Integrated with the Phoenix Micron Retinal Imaging Microscope and associated hardware	

**We have done our research so that you can do yours.**